

cardiofax^C

Electrocardiograph ECG-2150



Easy access to reliable diagnosis

Fighting Disease with Electronics

 **NIHON KOHDEN**

Easy access to reliable diagnosis

Reliable and accurate diagnosis

Cardiofax C uses the Nihon Kohden ECG analysis program ECAPS 12C. It provides simultaneous 12-lead ECG acquisition and analysis with approximately 200 findings with active noise suppression. ECAPS 12C can also find typical waveforms of Brugada-type electrocardiogram. Accurate data supports fast diagnosis.

Real time monitoring

12-lead waveforms can be confirmed at a glance on the angled display before recording. This improves the daily workflow by avoiding a repeated examination.





Paperless data management

The compact cardiofax C has built-in LAN capability. The ECG data can be reviewed and compared on a PC with Nihon Kohden's optional ECG data management software polaris.one. Paperless ECG management saves time as well as cost.

Highest mobility in daily workflow


The lightweight and compact design of cardiofax C and a long-time battery operation capability of more than 120 minutes allows highest mobility in the daily workflow.

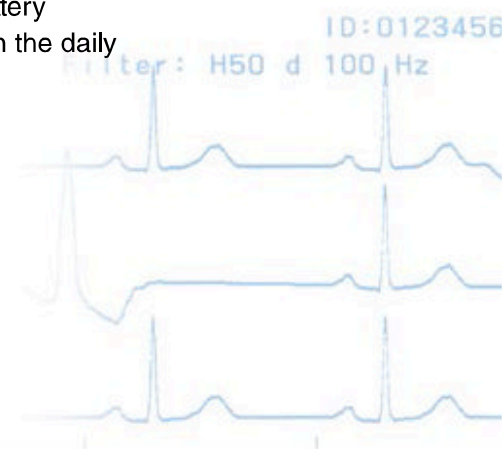

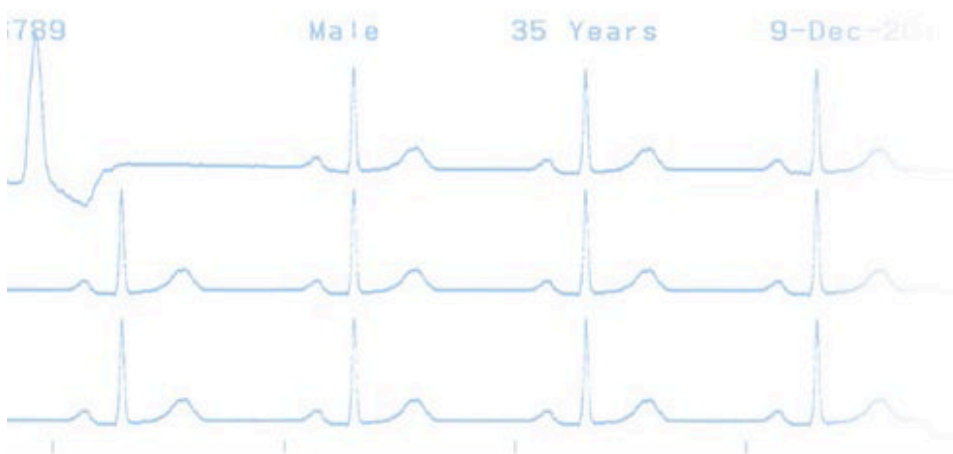


Recording Examples

12-lead ECG (3 traces × 4 groups) with analysis results

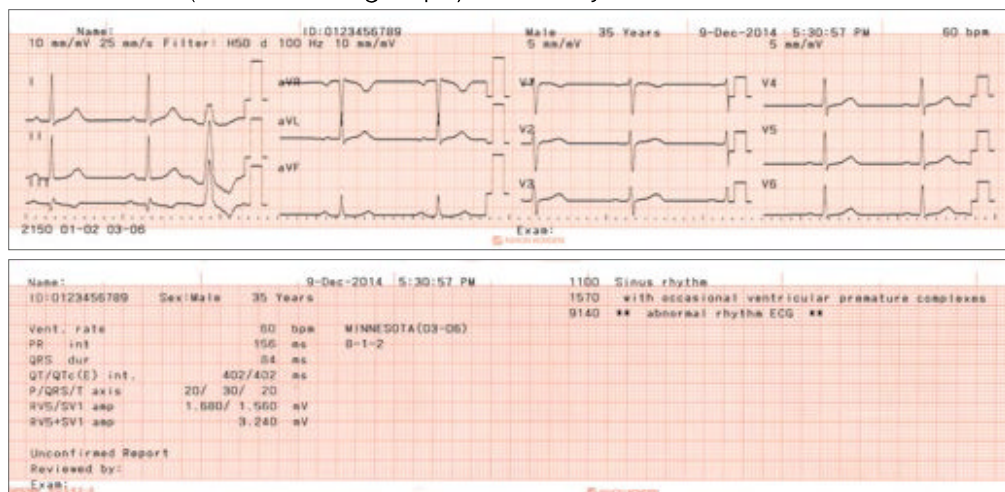

1 channel rhythm recording

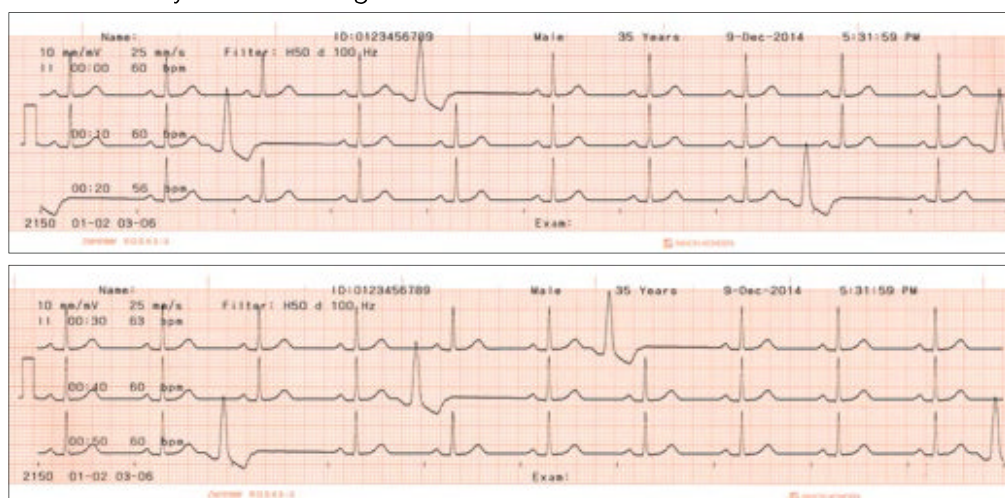

Specifications


Dimensions	279 W × 75 H × 210 D mm
Weight	1.7 kg (without battery and recording paper)
Battery operation time	120 minutes
Display	4.8-inch, 320 × 240 dots, monochrome LCD
Number of channels	3
Paper speed	25, 50 mm/s
Recording paper	63 mm width x 30 m, roll
Safety and performance standards	IEC 60601-1: 2012 IEC 60601-2-25: 2011 IEC 60601-1-2: 2007

Major Options

- ECG management software: polaris.one
- HL7 interface
- Barcode reader
- Magnetic card reader
- Cart, KD-105E
- Patient cable hanger, KH-100D



This brochure may be revised or replaced by Nihon Kohden at any time without notice.

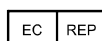

NIHON KOHDEN EUROPE GmbH
 Raiffeisenstrasse 10, 61191 Rosbach, Germany
 Phone: +49 (0) 6003 827-0, Fax: +49 (0) 6003 827-599
 Internet: www.nihonkohden.net, E-mail: info@nke.de



NIHON KOHDEN CORPORATION
 1-31-4 Nishiochiai, Shinjuku-ku, Tokyo 161-8560, Japan
 Phone: +81 (3) 59 96-80 36, Fax: +81 (3) 59 96-81 00
 Internet: www.nihonkohden.com