


Product Data Sheet:

I.V.-S86 - Non vented sets

S86-3W (dec., L-L)

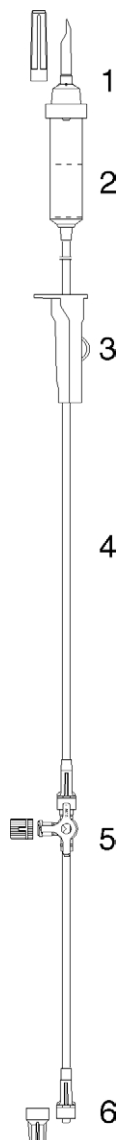

1. The special design of the spike enable easy penetration of container ports and offer a safe connection to the container. Non vented version. The protective cap maintains sterility of the set.
2. Standard drip tube (20 drops / 1 ml). The flexible, transparent non-reflective material with level ringmark enables easy priming and control of drip rate. With 15 micron fluid filter.
3. Precision roller clamp enables exact adjustment of the desired flow rate maintaining it with minimum variation. Furthermore clamp/hanging device for tubing/needle adapter and for spike protection after use of the set. (*) orange coloured roller clamp
4. Transparent flexible PVC tubing, i.d. 3.0 mm / o.d. 4.1 mm in the set (Regular bore).
5. Deconnectable 3-way stopcock with various flow way adjustments enables simultaneous administration medication with e.g. giving set or syringe.
6. Male Luer-Lock adapter guarantees a safe connection to the venipuncture device. FlowStop protective cap for male Luer-Lock adapters with integrated hydrophobic filter membrane, ensuring internal sterility prior to use. During priming the FlowStop cap remains attached to the Luer-Lock adapter allowing expulsion during priming of the set. The hydrophobic membrane prevents leakage. The risk of contamination is likewise reduced.

Ref.	Product Description	Packaging Units
37.2717	S86, 190 cm NoDEHP tube, deconnectable, 3-way stopcock, rotating Luer-Lock, FlowStop	80/320/2880
37.3160	S86, 210 cm NoDEHP tube, deconnectable, 3-way stopcock, Luer-Lock, FlowStop	100/3600